

Certificate of Analysis

SLAVIC LABS
SLAVIC LABS
STRAŽACKA 56
15-687 BIAŁYSTOK
Poland
michal@slaviclabs.com

Liquilabs s.r.o.
Inovační 122
252 41 Zlatníky-Hodkovice
Czechia
www.liquilabs.cz



[Verify Results Online](#)

Sample Identification

Sample Name Alarelin 5 mg
Batch Number
Date Published 2026-06-18 16:31

Results for LYO-0221

Peptides	Result	Unit	Uncertainty	Acceptable Range
Alarelin Assay Peptide Screening 0.1% TFA	4.93	mg	[± 0.02]	
Alarelin Purity Peptide Screening 0.1% TFA	> 99.8	%		
Alarelin Identification by Spectrum (FTIR) Peptide Screening 0.1% TFA	982		[± 5]	
Alarelin Identification by RT Peptide Screening 0.1% TFA	1.000		[± 0.005]	
Microbiology	Result	Unit	Uncertainty	Acceptable Range
Bacterial Endotoxin Chromgenic USP<85>/ Eur. Ph. 2.6.14. Bacterial Endotoxin Chromgenic Test	0.005	EU/mg	[± 0]	0 - 0.5

	Method Specification	
Determination of identity, content and purity of Alarelin		
<i>Document number</i> ALA_006_2026	<i>Superseded document</i> -	<i>Number of pages</i> 4

1. Content Assessment

1.1. Instrumentation

Module	Name	Serial Number
System Controller	Shimadzu CBM-40 Lite	L221226351398
Degassing Unit	Shimadzu DGU-403	NA
Pump	Shimadzu LC-40B XR	L22146350580
Autosampler	Shimadzu SIL-40C XR	L22216351622
Colum Thermostat	Shimadzu CTO-40S	L22236351602
PDA Detector	Shimadzu SPD-M40	L22276352808
SQ MS Detector	Shimadzu LCMS-2050	O12476200760

1.2. Chromatographic conditions

Chromatographic conditions	
Eluent A	0.05% TFA in Water (HPLC, Gradient Grade)
Eluent B	0.0425% TFA in Acetonitrile (HPLC, Gradient Grade)
Flow rate	0.9 mL/min
Program	Gradient elution
Injection volume	2 µL
Colum Temperature	55°C
Column	Waters XSelect CSH C18, 100x2.1mm 2.5µm
Detection wavelength	280nm


Gradient Program		
Time [min]	A [%]	B [%]
1.5	95	5
13	45	55
13.5	1	99
14.5	1	99
14.51	95	5
16	end	

1.3. Sample preparation

Whole amount of container was dissolved in 2mL of water (LCMS Grade). 100 µL of sample was transferred to HPLC vial and diluted by 900 µL water (LCMS Grade) and submitted for analysis.

1.4. Calibration curve

Calibration curve detail	
Quantitative method	External Standard
Calibration Type	Linear
Number of calibration points	5
Force through Zero	Enabled
Weighting Method	None

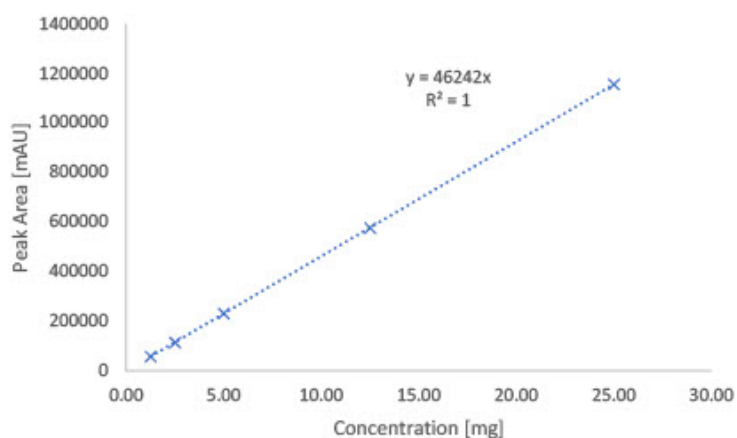

2. Purity assessment

2.1 Instrumentation

Module	Name	Serial Number
System Controller	Shimadzu CBM-40 Lite	L221226351398
Degassing Unit	Shimadzu DGU-403	NA
Pump	Shimadzu LC-40B XR	L22146350580
Autosampler	Shimadzu SIL-40C XR	L22216351622
Colum Thermostat	Shimadzu CTO-40S	L22236351602
PDA Detector	Shimadzu SPD-M40	L22276352808
SQ MS Detector	Shimadzu LCMS-2050	O12476200760

2.2 Chromatographic conditions

Chromatographic conditions	
Eluent A	0.05% TFA in Water (HPLC, Gradient Grade)
Eluent B	0.0425% TFA in Acetonitrile (HPLC, Gradient Grade)
Flow rate	0.9 mL/min
Program	Gradient elution
Injection volume	2 µL
Colum Temperature	55°C
Column	Waters XSelect CSH C18, 100x2.1mm 2.5µm
Detection wavelength	225nm

Gradient Program		
Time [min]	A [%]	B [%]
1.5	95	5
13	45	55
13.5	1	99
14.5	1	99
14.51	95	5
16	end	

2.3 Purity assesment

Purity of compound assesed by area normalization method, comparing area of each peak to sum of area of all peaks detected at wavelenght of 214 nm.

3. Identity Assessment

3.1 Instrumentation

Module	Name	Serial Number
System Controller	Shimadzu CBM-40 Lite	L221226351398
Degassing Unit	Shimadzu DGU-403	NA
Pump	Shimadzu LC-40B XR	L22146350580
Autosampler	Shimadzu SIL-40C XR	L22216351622
Colum Thermostat	Shimadzu CTO-40S	L22236351602
PDA Detector	Shimadzu SPD-M40	L22276352808
SQ MS Detector	Shimadzu LCMS-2050	O12476200760

3.2 Chromatographic conditions

Chromatographic conditions	
Eluent A	0.05% TFA in Water (HPLC, Gradient Grade)
Eluent B	0.0425% TFA in Acetonitrile (HPLC, Gradient Grade)
Flow rate	0.9 mL/min
Program	Gradient elution
Injection volume	2 µL
Colum Temperature	55°C
Column	Waters XSelect CSH C18, 100x2.1mm 2.5µm
Mass spectrometry	Scan: positive 280-2000 Da

Gradient Program		
Time [min]	A [%]	B [%]
1.5	95	5
13	45	55
13.5	1	99
14.5	1	99
14.51	95	5
16	end	

3.3 Molecular Ion Mass evaluation

Molecular ion mass was determined by deconvolution of multiply charged ESI-MS spectra to calculate the average neutral (zero-charge) molecular mass by equation:

$$M(\text{neutral}) = (z_i((mz_i) - H)) - ME$$

Where:

mz_i - Measured mass of charged particle

z_i - charge

H - proton mass (1.0076 Da)

ME – mass error

Analysis Report




Analysis of quantity and purity of active ingredient by UHPLC with UV detection

Sample Information


Injection Volume : 2
Data File : LYO-0221__006.lcd
Method File : Peptide screening_V7_Group B.lcm
Date Acquired : 6/12/2026 2:50:07 PM

Chromatogram

uAU

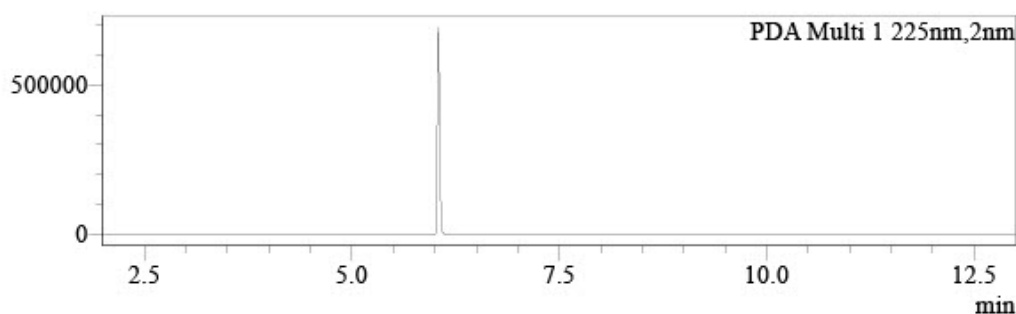

uAU

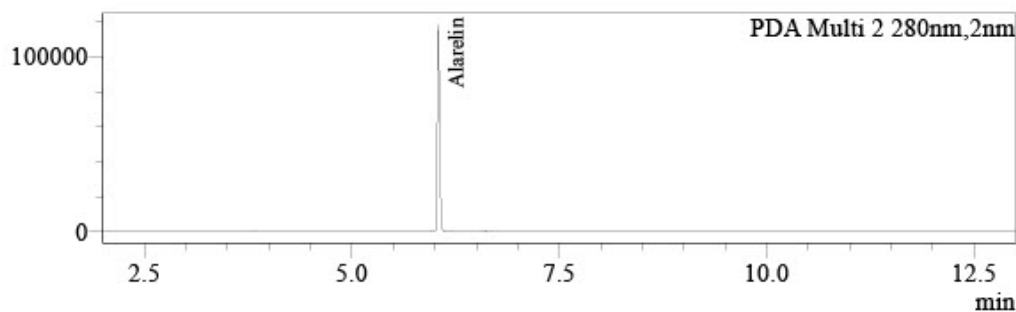

Peak Table

PDA Ch1 225nm

Name	Ret. Time	Area	Conc.	Unit	Area%
	4.253	469	0.000		0.034
	5.323	800	0.000		0.058
	5.963	115	0.000		0.008
	6.041	1370175	0.000		99.837
	6.228	47	0.000		0.003
	6.352	204	0.000		0.015
	8.631	87	0.000		0.006
	11.765	508	0.000		0.037
		1372405			100.000

Peak Table

PDA Ch2 280nm

Name	Ret. Time	Area	Conc.	Unit
Alarelin	6.041	227956	4.930	mg
		227956		

Analysis Report




LIQUILABS

Analytical services


Analysis of identity of active ingredient by UHPLC with mass spectrometric detection

MS Chromatogram

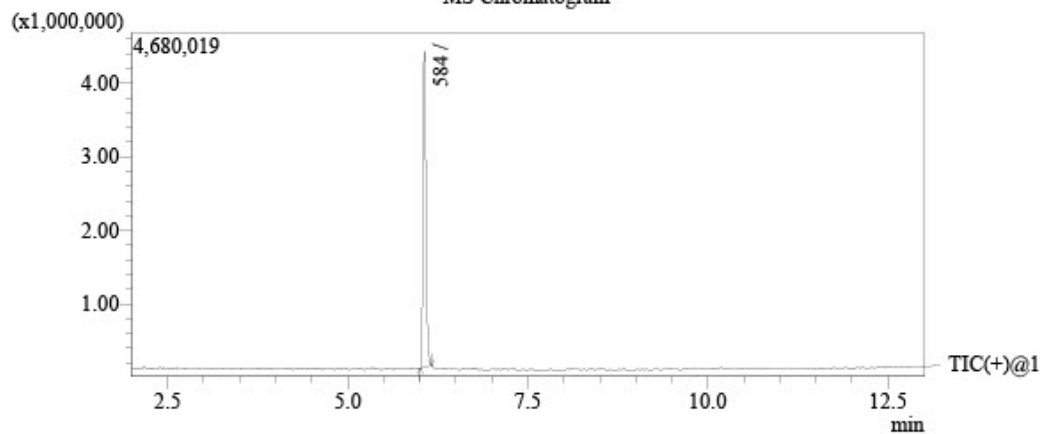


Library Search

Detected compound
Ret. Time: 6.008 Base Peak: 584
Scan Range: 280.00 - 2000.00 Polarity: Positive

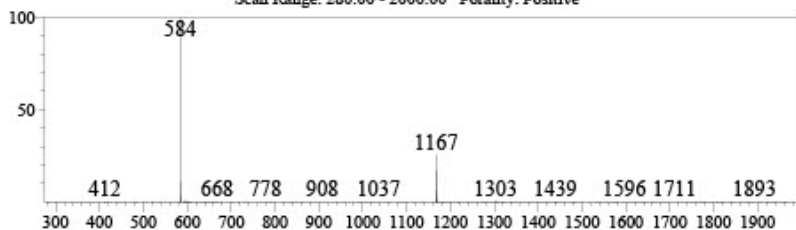

Reference

Library: Peptide screening lib
Formula: C₅₆H₇₈N₁₆O₁₂ CAS: 79561-22-1 Mol. Weight: 1167 Ret. Index: 0
Compound Name: Alarelin

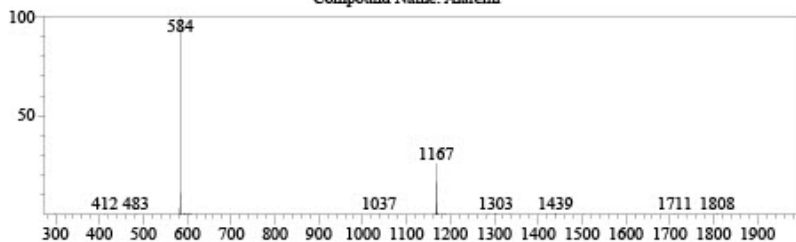

	Method Specification	
Determination of bacterial endotoxin content of lyophilized samples		
<i>Document number</i> ENDOTOX_0612_2026	<i>Superseded document</i> -	<i>Number of pages</i> 2

1. Chromgenic LAL Assay Determination of Bacterial Endotoxin content of sample

1.1. Instrumentation

- Pipette set 1-1000 µL
- Thermostatically controlled water bath
- UV VIS spectrometer (Shimadzu UV-1601)
- GenScript ToxinSensor Chromgenic LAL Endotoxin Assay kit

1.2. Chemicals

- LAL Reagent water (endotoxin free)
- Limulus Amoebocyte Lysate
- LAL Substrate
- Color Stabilizer #1
- Color Stabilizer #2
- Color Stabilizer #3
- 35% HCl (p.a.)

1.3. Sample preparation

1. Sample container was weighed prior to dissolution and measured weight was marked.
2. Sample was completely dissolved in its container by 2 mL of LAL Reagent water.
3. 100 µL of the sample was aliquoted for analysis.
4. After analysis container was emptied and dried.
5. Dry mass of container was measured and exact weight of dissolved content was determined as:

$$m_{dc} = m_{sample} - m_{container}$$

1.4. Toxin sensor Chromgenic LAL Endotoxin Assay kit preparation


Procedures regarding preparation of reaction solutions possible to find in:

https://www.genscript.com/site2/document/5292_20080806231827.PDF

1.5. Measurement procedure

	Standards	Samples	Blank
Standards (mL)	0.1	-	-
Samples (mL)	-	0.1	-
LAL Reagent Water (mL)	-	-	0.1
LAL Solution (mL)	0.1	0.1	0.1
Mix well and incubate at 37°C for 27 min			
Substrate solution (mL)	0.1	0.1	0.1
Mix well and incubate at 37°C for 6 min			
Color Stabilizer #1 solution	0.5	0.5	0.5
Color Stabilizer #2 solution	0.5	0.5	0.5
Color Stabilizer #3 solution	0.5	0.5	0.5
Mix well and read the absorbance at 545nm			

1.6. Calibration curve

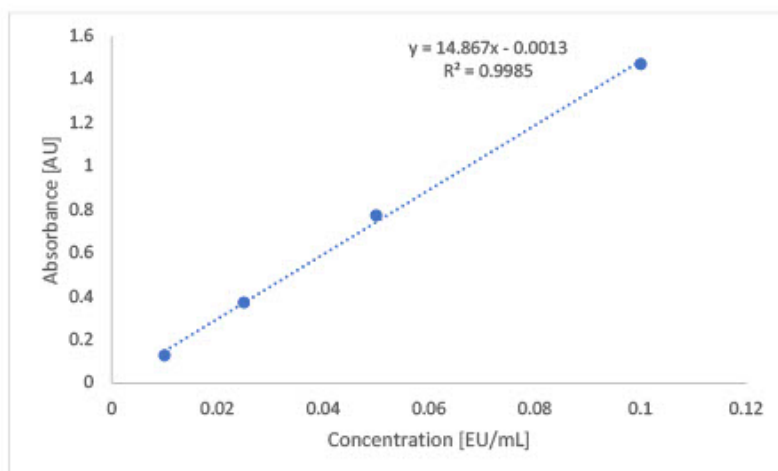

1.7. Calculation of endotoxin content

Endotoxin content of the sample was calculated from the calibration curve as:

$$Endotox[EU/mg] = \frac{\left(\frac{ABS_{sample}}{S_{calib}} \right) * 20}{m_{sample}}$$

ABS_{sample} = Measured absorbance of sample

S_{calib} = Slope of calibration curve

m_{sample} = real measured mass of sample

20 = dilution factor of measured sample

Responsibles



Mr. Ján Galbavý
CEO

Analysis results relate only to the samples tested.

This document shall not be reproduced except in full, without the written approval of Liquilabs s.r.o.